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Thank you for purchasing our product.
Please read all information before using this product.
Please keep this user manual for future reference.

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WARNINGS

- 1) Keep out of reach of children. Do not allow children to use without adult supervision.
- Noncompliance with this instruction may in children swallowing the toothbrush.
- If swallowed by children, consult a physician immediately.
- 2) This product is for personal use only and is not intended professional use in dental office or institution.
- 3) Consult your dental professional before using this product.
Also please do not use if you are
- hypersensitive in your mouth
- during treatment and after periodontal operation
- having treatment or postoperative in your mouth
- bleeding in your mouth
- 4) Do not use this product with head tilted backwards or if the brush comes off, there is a possibility of suffocation.
- 5) Do not modify this product.
Noncompliance with this instruction may result in injury or damage.
- 6) Persons with pacemakers or other implanted devices should consult with a physician about any adverse affects on their pacemaker that may result from using this product.
When operated, this product may affect pacemakers.
- 7) Please use "AAA alkaline battery".
- 8) Ensure that young children do not play with the product.
- 9) Do not brush on the tooth and gum with strong pressure.
And do not put plastic part, except bristles, on your teeth and gum.
Also do not brush the same spot for a long time.
A tooth and the gums can be injured and cause the breakage.
- 10) Stop using when you feel uncomfortable.
Discontinue use if this product appears damaged in any way (brush head or main body).
- 11) Avoid using a brush head with crushed or bent bristles.
Damaged brush may break off while brushing.
- 12) This product is only for personal teeth and gum brushing use only.
- 13) Keep the power switch off when not in use.
- 14) Make sure that the brush is attached securely.
- 15) While it is inside your mouth, please hold the brush properly.
- 16) Keep the main body away from water.
- 17) Do not use this product in water or shower.
- 18) Do not turn the power switch ON before putting this product into your mouth.
Turn the power switch OFF before taking out the product from your mouth.
- 19) Do not use toothpaste containing peroxide, baking soda, or whitening agents.
- 20) Do not use the product for other than intended use.

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HOW TO CARE


Soak the brush part in warm water to get rid of toothpaste residue between bristles.

Do not immerse the main body and the battery cover in water.

*In case of getting dirty on the main body, wipe surface with a damp cloth immediately.

*This product should never be placed in the dishwasher.

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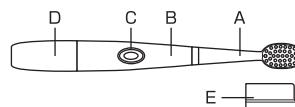
Hold the main body and take the battery cover off as shown below.

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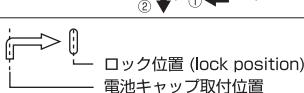
Insert a new battery with the(-)side facing up.

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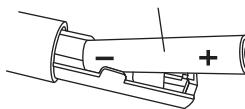
*Make sure the polarity of the battery is correctly aligned.

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Hold the main body and put the battery cover on as shown below.

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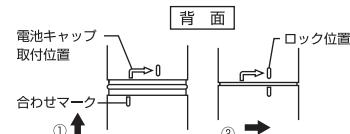
After putting the brush in the mouth, press the power switch.

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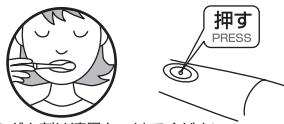※Put an appropriate amount of toothpaste on the brush.

※Press the power switch until vibrations begin. (Push the power switch gently during operation to pause)

※Do not press the power switch with a claw. Rubber of the power switch may get stuck.

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Press the power switch until vibrations stop after you finish brushing.

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※Be sure the power is off before removing the toothbrush from your mouth, otherwise water and toothbrush may splatter.

※Press the power switch until vibrations stop. (Pressing the switch gently pauses but releases your finger it will start moving again.)

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Rinse the brush head and brush with water after each use.
Do not immerse the main body in water.

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
※Let the brush head air-dry before storing.

※Make sure the battery box is secure to prevent water from getting inside the main body.

ブラシの交換 How to replace brush

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Hold the main body and remove the brush.

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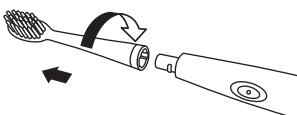
※Turn the brush counterclockwise, then pull it off the main body.

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Place the new brush head on the main body as shown below.

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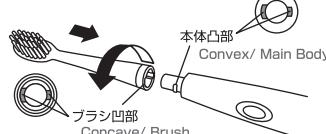※Turn the brush head clockwise until the brush is securely attached on the main body.

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The brush heads should be replaced every 3 months to achieve optimal results.(Sold separately).

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
歯の面に対してブラシを直角にあて、1~2歯ずつヨコ方向へ動かしてみがく。(5往復)
Place the brush perpendicularly against the tooth surface and move the brush toward one to two teeth of the side. (5 round trips)

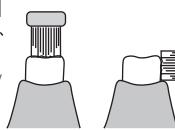歯と歯ぐきの境目 (between teeth and gums)

ブラシの毛先が歯ぐきの溝に入るよう斜めに傾け(約45度)1~2歯ずつヨコ方向へ動かしてみがく。(5往復)
Tilt diagonally so that the tip of the brush enters the groove of the gum. (About 45°) Move the brush toward one to two teeth of the side. (5 round trips)

歯の裏や歯並びの悪い箇所 (place on the back of teeth and bad placement of teeth)

歯ブラシをタテにあて、歯の面に沿って1本ずつみがく。
Place the toothbrush vertically and brush one by one along the tooth surface.

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